

SoftCompress

Compression bandages

When compression bandaging is used, padding underneath prevents pressure points and skin damage. It can also improve quality of life¹. In addition, the SoftCompress cushioning material can significantly increase the success of decongestive therapy. This is thanks to the vertical foam channels that form a unique surface structure to support lymphatic drainage. Bandaging aids and additional products feature a particularly skin-friendly cotton surface, as opposed to frequently used synthetic materials. This provides the patient with pleasant wearing comfort. In addition, SoftCompress bandaging aids make it easy for patients to use the bandages themselves, therefore ensuring the long-term success of treatment.

Features

- + core made of high-quality foam
 - optimum transfer of pressure
 - perfect pressure distribution
- + vertically-stitched foam channels
 - increased lymph drainage
 - effective for counteracting fibrosis
- + exceptionally skin-friendly, thanks to the cotton surface
 - breathable
 - smooth
- + constriction in the bends of the joints is prevented
- + protects against bandage edges pinching
- + saves time when bandaging
- + easier self-bandaging
- + washable up to 60 °C
- + suitable for tumble-drying at a low temperature
- + disinfection and treatment possible – therefore saves money²
- + **SoftCompress bandages can be worn under shortstretch bandages or in combination with adjustable compression wraps**

Indications

- lymphoedema
- lipoedema
- phleboedema
- post-traumatic oedema
- postoperative oedema
- fibrosis

Material composition

- Top cover: 100 % cotton
- Foam: 100 % polyurethane

Contraindications

See page 37

Colour



White



With cotton



Washable
60 °C



Suitable for
tumble-drying



Foam channels

Vertically-stitched foam channels promote lymphatic drainage and help to loosen fibrosis.



Effectiveness

Visible effects of SoftCompress bandages in cases of secondary leg lymphoedema.

¹Protz et al. "Kompressionsbandagierungen mit und ohne Unterpolsterung im Vergleich" (Comparing compression bandaging with and without cushioning). Vasomed 2019; 31(3) 132-133.

²The products have been successfully tested in accordance with DIN ISO TS 15883-5, SOD 19-002, EN ISO 17664, SOD 19-01 and the Robert Koch Institute "Hygiene Requirements for the Reprocessing of Medical Devices" guideline.

Bandages

In phase 1 of CDT, the comfortable SoftCompress bandages can be worn under short-stretch bandages or, in phase 2, under medical compression stockings or ACS Light. They provide additional compression, cushioning and an ideal pressure distribution. The bandages enable bandaging to be carried out more quickly and effectively. They also make it easier for patients to apply the bandages themselves. They are available in uni size, standard size or custom-made. Custom-made bandages can also be used for treatment at night.

For CDT phase 1



**SoftCompress
Bandage Lower Leg**
in uni size
Art. 6831



**SoftCompress
Bandage Lower Leg and Thigh**
in uni size
Art. 6821



**SoftCompress
Bandage Arm**
in uni size
Art. 6823

Uni size

Uni size bandages are open on the back and are secured in place with the help of touch fasteners or an elastic fabric sleeve. This enables the bandaging aid to accommodate large fluctuations in circumference during therapy.

Standard and uni sizes can be trimmed at the markings provided.

Caution: Before trimming between the cross-seams, the product has to be washed once and stretched lengthwise while it's still wet, before drying.



For CDT phase 2



**SoftCompress
Bandage Lower Leg**
in standard sizes or custom-made
Art. 6830









**SoftCompress
Bandage Lower Leg and Thigh**
custom-made
Art. 6820



**SoftCompress
Bandage Arm**
in standard sizes or custom-made
Art. 6822


Product range

Art.	Model	Standard sizes	Custom-made	Uni size
 6831	Bandage lower leg (for CDT phase 1) Incl. 2 malleo pads	–	•	•
 6830	Bandage lower leg (for CDT phase 2)	1 – 5	•	–
 6821	Bandage lower leg and thigh (for CDT phase 1) Incl. bandage lower leg, uni size 2 malleo pads 1 knee pad 1 compression roll (10 cm x 4 m) 1 compression roll (12 cm x 4 m)	–	–	•
 6820	Bandage lower leg and thigh (for CDT phase 2) Incl. 2 malleo pads 1 knee pad	–	•	–
 6823	Bandage arm (for CDT phase 1) Incl. 1 hand pad 1 elbow pad	–	–	•
 6822	Bandage arm (for CDT phase 2) Incl. 1 hand pad 1 elbow pad	1 – 6	•	–

CDT= Complex decongestive therapy


Sizing charts

Juzo SoftCompress Bandage Lower Leg

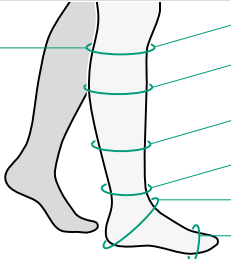



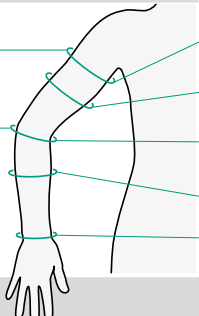

Art. 6830 Length in cm	Standard sizes - Circumferences in cm				
	1 / XS	2 / S	3 / M	4 / L	5 / XL
 PD 38.0 – 46.0	cD 27.0 – 29.5	29.5 – 32.5	32.5 – 35.5	35.5 – 39.0	38.5 – 42.0
	cC 29.0 – 31.5	31.5 – 34.5	34.5 – 37.5	37.5 – 41.0	40.5 – 44.5
	cB' 23.5 – 25.5	25.5 – 28.0	27.5 – 30.0	30.0 – 33.0	32.5 – 35.5
	cB 18.5 – 20.0	20.0 – 22.0	22.0 – 24.0	24.0 – 26.0	26.0 – 28.5
	cY 28.0 – 30.5	29.0 – 31.5	31.0 – 34.0	33.0 – 36.0	35.0 – 38.5
	cA 18.5 – 20.0	20.0 – 22.0	22.0 – 24.0	24.0 – 26.0	26.0 – 28.5

Important:
Standard sizes
are determined by
skin dimensions. For
custom-made
products, please
measure the item
"on tension".

Juzo SoftCompress Bandage Arm

Art. 6822 Lengths in cm	Standard sizes - Circumferences in cm					
	1 / XS	2 / S	3 / M	4 / L	5 / XL	6 / XXL
 PG 46.0 – 51.0	cG 30.0 – 33.0	32.0 – 35.0	34.0 – 37.0	36.0 – 39.5	38.0 – 41.5	40.0 – 44.0
	cF 28.0 – 30.5	30.0 – 33.0	32.0 – 35.0	34.0 – 37.0	36.0 – 39.5	38.0 – 41.5
PE 21.0 – 25.0	cE 25.0 – 27.5	27.0 – 29.5	29.0 – 31.5	30.5 – 33.5	32.0 – 35.0	33.5 – 36.5
	cD 23.0 – 25.0	24.5 – 27.0	26.0 – 28.5	27.5 – 30.0	29.0 – 31.5	30.5 – 33.5
	cC 14.5 – 16.0	16.0 – 17.5	17.5 – 19.0	19.0 – 20.5	20.5 – 22.0	22.0 – 24.0

Product range

Art.	Model	Uni size
	6800 KL LymphAid Leg, small (for CDT phase 1)	
	1 lower leg bandage, uni size 20 Rapidelast V – gauze bandages 4 cm x 4 m 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 8 cm x 5 m 2 short-stretch bandages 10 cm x 5 m 4 short-stretch bandages 12 cm x 5 m 2 malleo pads	.
	6800 GR LymphAid Leg, large (for CDT phase 1)	
	1 lower leg bandage, uni size 1 compression roll 10 cm x 4 m 1 compression roll 12 cm x 4 m 20 Rapidelast V – gauze bandages 4 cm x 4 m 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 8 cm x 5 m 4 short-stretch bandages 10 cm x 5 m 6 short-stretch bandages 12 cm x 5 m 2 malleo pads 1 knee pad	.
	6801 LymphAid Arm (for CDT phase 1)	
	1 arm bandaging aid, uni size 20 Rapidelast V – gauze bandages 6 cm x 4 m 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 10 cm x 5 m 3 short-stretch bandages 12 cm x 5 m 1 hand pad 1 elbow pad	.
	6810 KL BandageAid Leg, small (for CDT phase 1)	
	1 compression roll 10 cm x 4 m 20 Rapidelast V – gauze bandages 4 cm x 4 m 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 3 short-stretch bandages 10 cm x 5 m 2 malleo pads	.
	6810 GR BandageAid Leg, large (for CDT phase 1)	
	1 compression roll 10 cm x 4 m 2 compression rolls 12 cm x 4 m 20 Rapidelast V – gauze bandages 4 cm x 4 m 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 3 short-stretch bandages 10 cm x 5 m 4 short-stretch bandages 12 cm x 5 m 2 malleo pads 1 knee pad	.
	6811 BandageAid Arm (for CDT phase 1)	
	1 compression roll 10 cm x 4 m 20 Rapidelast V – gauze bandages 6 cm x 4 m 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 5 short-stretch bandages 10 cm x 5 m 1 hand pad 1 elbow pad	.



Application videos:

LymphAid, lower leg



LymphAid, leg



BandageAid, arm and hand



LymphAid, lower leg



BandageAid, leg



BandageAid, arm and hand



CDT= Complex decongestive therapy

The gauze bandages and short-stretch bandages used in LymphAid and BandageAid are quality products manufactured by Karl Otto Braun GmbH & Co. KG (KOB). They are also available separately.

Compression roll and compression sheet

SoftCompress bandaging materials are used for lymphological and phlebological compression bandaging. Thanks to their soft surface, they are very smooth, skin-friendly and breathable.



SoftCompress Compression Roll

Art. 6835

The thin and smooth SoftCompress Compression Rolls exert compression when applied correctly. The compression rolls evenly distribute the pressure of the short-stretch bandage applied on top, thus preventing constriction.



SoftCompress Compression Sheet

Art. 6836

The SoftCompress Compression sheets can be cut to size and can therefore be used in a number of different body regions. We recommend stitching down the edges.


Product range

Art.	Model	Standard sizes
6835	Compression Roll (for CDT phase 1)	
	5 mm, 8 cm x 4 m 5 mm, 10 cm x 4 m 5 mm, 12 cm x 4 m 5 mm, 16 cm x 4 m	• • • •
6836	Compression Sheet (for CDT phase 1)	
	65 x 65 cm 125 x 100 cm	• •


CDT= Complex decongestive therapy

Sizing charts

Juzo SoftCompress Compression Roll

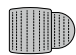Art. 6835	Standard sizes	Measurements in cm (length x breadth)
Thickness in mm: 5.0 	1	400.0 x 8.0
	2	400.0 x 10.0
	3	400.0 x 12.0
	4	400.0 x 16.0


Juzo SoftCompress Compression Sheet

Art. 6836	Standard sizes	Measurements in cm (breadth x height)
	1	65.0 x 65.0
	2	125.0 x 100.0

Additional products

SoftCompress


SoftCompress pads can be worn underneath compression bandages or compression garments. They are especially suited to providing a targeted increase in pressure in localised oedemas.

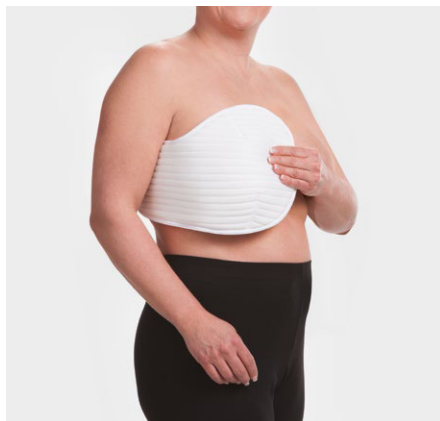

SoftCompress Breast Pad

Art. 6837

The SoftCompress Breast Pad serves as a comfortable pad for a compression bra or thorax compression vest and can be used for breast or chest wall oedemas, for example.

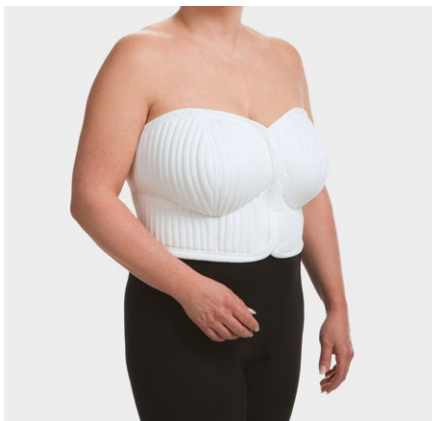

SoftCompress Breast Bandage

Art. 6838

The SoftCompress Breast Bandage combines the positive characteristics of the SoftCompress material with those of a compression bra. The bandage has an adjustable chest circumference. The breast bandage can be used for breast or chest wall oedema, for example.



SoftCompress Genital Pad for women

Art. 6839

The SoftCompress Genital Pad for women is inserted into compression tights. It provides a targeted increase in pressure, cushioning and an ideal pressure distribution when treating genital oedemas.



SoftCompress Genital Pad for men

Art. 6840

The SoftCompress Genital Pad for men (reversible) is inserted into compression tights. The wide channel structure helps to increase pressure and lymph drainage, while the narrow channels are useful for reducing hardening caused by fibrosis. It also provides a targeted increase in pressure and an ideal pressure distribution when treating genital oedemas.



SoftCompress Knee Pad

Art. 6845

This SoftCompress pad provides cushioning for the back of the knee joint.



SoftCompress Elbow Pad

Art. 6842

This SoftCompress pad provides cushioning where needed. It is especially useful for the crook of the elbow.



SoftCompress Malleo Pad

Art. 6843

SoftCompress Malleo Pads cushion the Bisgaard (perimalleolar) region of the ankle. This ensures a targeted increase in pressure and an ideal pressure distribution. The pads come in pairs.



SoftCompress Hand Pad

Art. 6844

SoftCompress Hand Pads are specially designed for the back of the hand and feature two finger openings. They ensure a targeted increase in pressure and an ideal pressure distribution.



SoftCompress Short-stretch Bandage

Art. 6846

The short-stretch bandages are made from 100 % cotton, creating compression bandaging with low resting pressure and high working pressure. They can be washed up to 95 °C and sterilised with steam (121 °C).



SoftCompress Rapidelast V gauze bandages

Art. 6847/6848

The Rapidelast V gauze bandages are made from 60 % polyamide and 40 % viscose. They are used to secure compression and cushioning material – especially on sensitive skin or joints – without restricting the patient's freedom of movement. Bandages with Rapidelast V are non-slip. The high elasticity of the fabric makes application quick and easy.



SoftCompress Tubular Bandages

Art. 6849

The tubular bandage is made from 100 % cotton. It is particularly stretchy and gentle on the skin. It protects the skin without constricting it and can also be used to secure pads.

Product range

Art.	Model	Standard sizes	Custom-made
6837	Breast pad (for CDT phases 1 and 2)		
	Cup A – B Cup C – D Cup E – F	1 – 3	•
6838	Breast bandage (for CDT phases 1 and 2)		
	Cup A Cup B Cup C Cup D Cup E Cup F	1 – 6	•
6839	Genital pad for women (for CDT phases 1 and 2)		
	20 x 20 cm 30 x 28 cm 40 x 32 cm	1 – 3	•
6840	Genital pad for men (for CDT phases 1 and 2)		
	20 x 28 cm 23 x 30 cm 26 x 32 cm	1 – 3	•
6842	Elbow pad (for CDT phases 1 and 2)	Uni size	
	10 x 13 cm (custom-made max. 13 x 16 cm)	•	•
6843	Malleo pads (for CDT phases 1 and 2)	Uni size	
	6.5 x 9 cm delivered in pairs	•	•
6844	Hand pad (for CDT phases 1 and 2)	Uni size	
	8.5 x 24 cm	•	•
6845	Knee pad (for CDT phases 1 and 2)	Uni size	
	13 x 17 cm (custom-made max. 19 x 23 cm)	•	•


Art.	Model	Each	Box
6846	Short-stretch bandage (for CDT phase 1)		10 items
	6 cm x 5 m 8 cm x 5 m 10 cm x 5 m 12 cm x 5 m	• • • •	• • • •
	Rapidelast V gauze bandage (for CDT phase 1)		20 items
6847	6 cm x 4 m (for the fingers)	–	•
6848	4 cm x 4 m (for the toes)	–	•
6849	Tubular bandage (for CDT phase 1)		
	6 cm x 20 m 8 cm x 20 m 10 cm x 20 m	• • •	– – –

Sizing charts


Juzo SoftCompress Breast Pad

Art. 6837	Standard sizes	Cup
	1	A – B
	2	C – D
	3	E – F


Juzo SoftCompress Breast Bandage

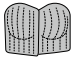Art. 6838	Standard sizes	Cup	Breast circumference in cm	Underbreast circumference in cm
	1	A	90.0 – 106.0	84.0 – 100.0
	2	B	90.0 – 106.0	84.0 – 100.0
	3	C	95.0 – 112.0	85.0 – 102.0
	4	D	95.0 – 112.0	85.0 – 102.0
	5	E	100.0 – 116.0	88.0 – 104.0
	6	F	100.0 – 116.0	88.0 – 104.0

Juzo SoftCompress Genital Pad Women

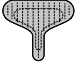Art. 6839	Standard sizes	Measurements in cm (breadth x height)
	1	20.0 x 20.0
	2	30.0 x 28.0
	3	40.0 x 32.0

Juzo SoftCompress Genital Pad Men

Art. 6840	Standard sizes	Measurements in cm (breadth x height)
	1	20.0 x 28.0
	2	23.0 x 30.0
	3	26.0 x 32.0